

HAEGARDA® Reconstitution Guide

Haegarda and Mix2Vial®
work together to help you
with self-administration

1 Clean Stoppers

Remove the flip caps from both vials (HAEGARDA and diluent). Wipe rubber stoppers with an antiseptic wipe and allow to dry.



3 Prepare Diluent Vial

Place the diluent vial on an even, clean surface and grip it tightly. Take the Mix2Vial® together with its package and push the spike of the blue adapter end straight down through the diluent vial stopper.



5 Transfer Diluent into Haegarda Vial

Place the HAEGARDA vial on a firm, even surface. Invert the diluent vial with the Mix2Vial® set attached, and push the spike of the transparent end of the Mix2Vial® straight through the HAEGARDA vial stopper. The diluent will automatically flow into the HAEGARDA vial.



2 Open the Mix2Vial®

Open the Mix2Vial® package by peeling off the lid, leaving the Mix2Vial® inside the package to maintain sterility. Do not remove the Mix2Vial® from the blister package!



4 Remove the Mix2Vial® Packaging

While holding the diluent vial, carefully remove the package from the Mix2Vial® set by pulling it vertically upwards. Make sure to pull away only the package and not the Mix2Vial® set.

BEFORE YOU BEGIN

ASSEMBLE THE SUPPLIES PROVIDED IN THE PRODUCT PACKAGE ON A FLAT SURFACE (INFUSION MAT):

- HAEGARDA and diluent vials (must be at room temperature)
- Mix2Vial® needleless transfer device
- Sterile syringe
- Subcutaneous infusion set or hypodermic needle
- Alcohol or disinfectant wipes (not included in product package)

Thoroughly clean your working surface using alcohol or disinfectant wipes, and thoroughly wash and dry your hands.

Use aseptic technique when preparing and withdrawing HAEGARDA.



6 Dissolve HAEGARDA

With the diluent and HAEGARDA vials still attached to the Mix2Vial® set, gently swirl the HAEGARDA vial to ensure that the powder is fully dissolved (generally within 5 minutes but may take as long as 10 minutes). Do not shake the vial.



8 Load the Syringe

Draw air into an empty, sterile syringe. With the HAEGARDA vial upright, screw the syringe clockwise into the attached Mix2Vial® set. Inject air into the HAEGARDA vial.



7 Unscrew Empty Diluent Vial


Grasp the transparent HAEGARDA side of the Mix2Vial® set with one hand, and grasp the blue diluent side with the other hand. Unscrew the set carefully counter-clockwise into two pieces. Discard the diluent vial with the blue Mix2Vial® adapter attached.



9 Draw Solution into Syringe

Keeping the syringe plunger pressed, turn the system upside down and draw the solution into the syringe by pulling the plunger back slowly.

For more information,
please visit
patients.cslbehring.ca

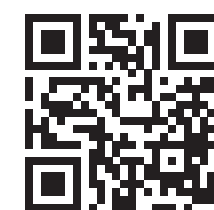

10 Prepare Administration Set

Once the solution has been transferred into the syringe, firmly hold the barrel of the syringe (keeping the plunger facing down) and unscrew it counter-clockwise from the Mix2Vial® set. The reconstituted solution should be colourless, clear and free of visible particles. Do not use if the solution contains particulates or is discoloured. As instructed by a healthcare provider, attach a hypodermic needle or subcutaneous infusion set. Prime the needle or tubing as required and instructed.



Record the product lot number, date and time of infusion in your treatment log book after every use of HAEGARDA.

For complete information and administration steps, please refer to the Product Monograph available at <https://www.cslbehring.ca/products/product-list>.

Mix2Vial® is a registered trademark of West Pharma. Services IL, Ltd., a subsidiary of West Pharmaceutical Services, Inc. HAEGARDA® is a registered trademark of CSL Behring GmbH. HAEGARDA® is manufactured by CSL Behring GmbH and distributed by CSL Behring Canada Inc. © 2024 CSL Behring. All rights reserved.

HGD2024060179EN